
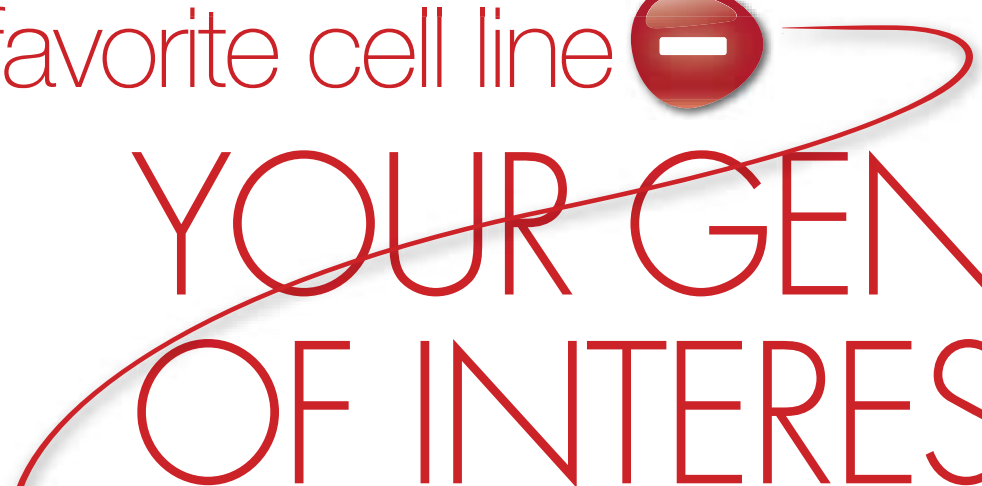


Your favorite cell line  
**YOUR GENE
OF INTEREST**



ADD.



DELETE.



SWAP.

EDIT FOR THE BETTER.


The next scientific breakthrough is here, and it's a knockout!

CompoZr™ ZFN technology, a breakthrough that enables simple, efficient and target specific genomic editing, exclusively from Sigma Life Science, is now available pre-designed and validated for knockout of your gene of interest. CompoZr Knockout ZFNs have been designed and validated for many human gene targets. Visit sigma.com/zfn to permanently knockout your gene of interest today.

Now it's all possible.



Sigma and Sigma-Aldrich are registered trademarks belonging to Sigma-Aldrich Co. and its affiliate Sigma-Aldrich Biotechnology, L.P.
CompoZr is a trademark of Sigma-Aldrich Co. and its affiliate Sigma-Aldrich Biotechnology, L.P.

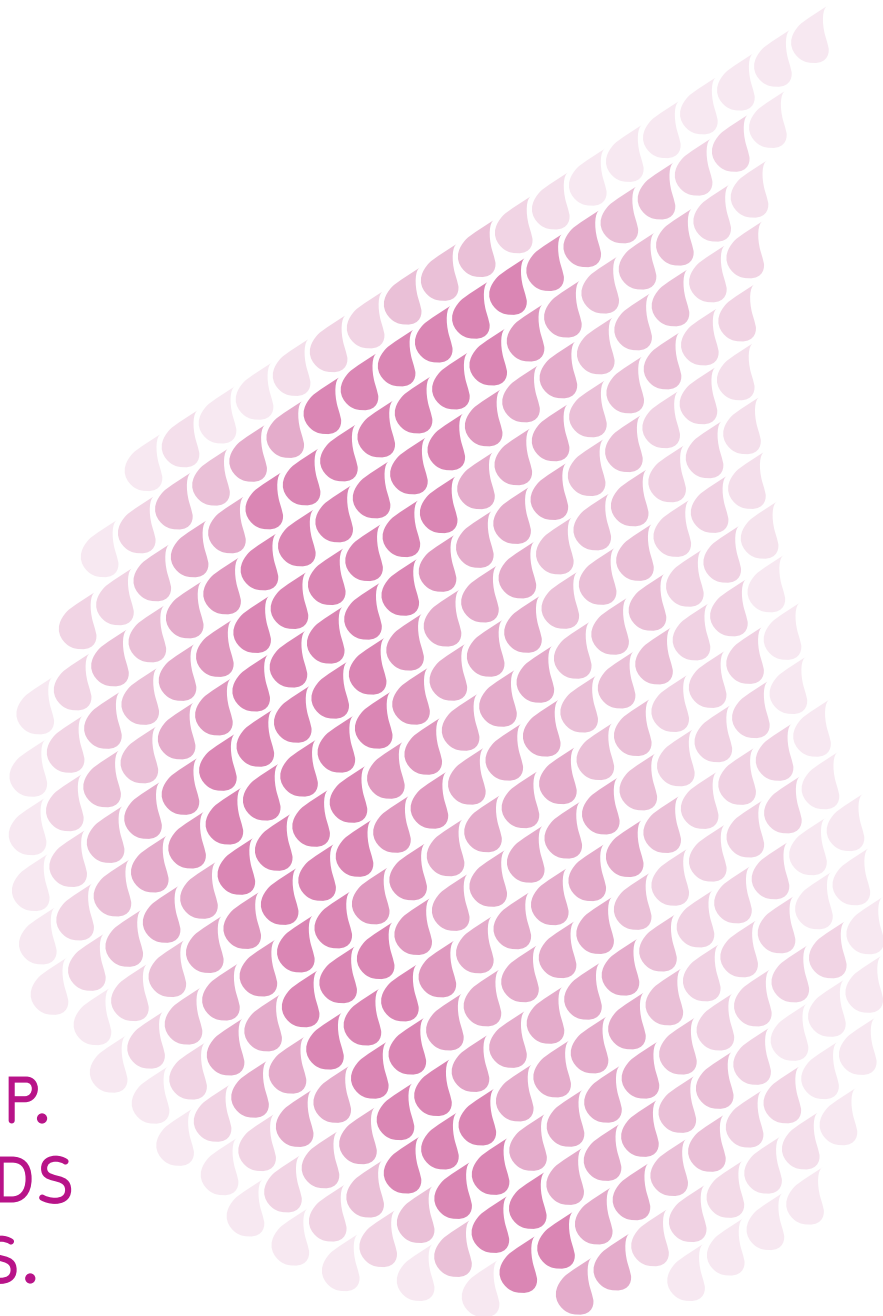

SOME SEE A DROP.
WE SEE HUNDREDS
OF EXPERIMENTS.

Choose Amicon® Ultra 0.5 mL for sample concentration so good,
even one drop will be enough to get all the downstream results you envision.

It's an ingenious design—a vertical filter that concentrates your smallest sample in as little as 10 minutes. Simply flip the filter upside down and spin it out to recover >90% with reliability. When downstream quality matters, make sure your upstream tools are the best. You'll see the difference.



ADVANCING LIFE SCIENCE TOGETHER™
Research. Development. Production.



Choose the perfect filter at
www.millipore.com/onedropB

Millipore and Amicon are registered trademarks of Millipore Corporation.
Advancing Life Science Together and the Millipore logo are trademarks of Millipore Corporation.
©2009 Millipore Corporation. All rights reserved.

Äsk

process developers what they
want in a protein purification
system and what do you get?



imagination at work

Pure excellence, pure efficiency

QIAxcel and QIAcube



- Automate your spin column preps on the QIAcube
- Eliminate tedious DNA fragment and RNA analysis with the QIAxcel
- Free up your time with automated purification and analysis

Visit www.qiagen.com/goto/PureExcellence for more information!






www.roche-applied-science.com

LightCycler® 480 Real-Time PCR System

The New Horizon in Real-Time PCR

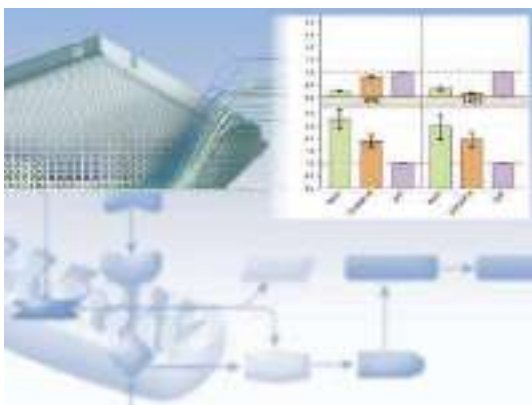

RealTime ready

Use these function-tested, custom qPCR assays and panels to add convenience and confidence to your gene expression analysis.

For life science research only. Not for use in diagnostic procedures.

This LightCycler® 480 Real-Time PCR System is a real-time thermal cycler licensed for research use only under U.S. Patent No. 6,814,934 and corresponding claims in its non-U.S. counterparts, and under one or more of U.S. Patents Nos. 5,475,510, 5,802,766, 6,703,296, 7,238,517, 7,594,241, 7,557,577 or corresponding claims in their non-U.S. counterparts, owned by Applied Biosystems, LLC. No rights are conveyed expressly, by implication or estoppel to any other patent claims. The LightCycler 480 Real-Time PCR System is not cleared for in vitro diagnostic use in the U.S. For further information on purchasing additional rights, contact the Director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California, 94404, USA. No right is conveyed expressly, by implication, or by estoppel under any other patent claim, such as claims to apparatus, reagents, kits, or methods such as 5' nuclease methods. License rights to practice PCR methods under the foregoing Roche patents for research and other non-in vitro diagnostic applications may be purchased from Applied Biosystems or may be obtained by purchasing licensed reagents from Roche, Applied Biosystems, or other Authorized third party. For information on purchasing licenses for research and other non-in vitro diagnostic applications, contact the Director of Licensing at Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

A license to practice PCR methods with other real-time detection under patents of F. Hoffmann-La Roche Ltd. and other Roche Molecular Systems („Roche“) for in vitro diagnostic applications may be purchased from Roche or may be obtained by purchasing in vitro diagnostic reagents from Roche or any other Authorized third party. The LightCycler 480 Real-Time PCR System is not cleared for in vitro diagnostic use in the U.S.

Practice of the patented 5' Nuclease Process requires a license from Applied Biosystems. The purchase of this product includes an immunity from suit under patents specified in the product insert to use only the amount purchased for the purchaser's own internal research when used with the separate purchase of Licensed Probe. No other patent rights are conveyed expressly, by implication, or by estoppel. Further information on purchasing licenses may be obtained from the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA.

Experimental data in this article are unpublished, courtesy by U. Tschulena, Ph.D., DKFZ Heidelberg, Germany.

© 2009 Roche Diagnostics GmbH. All rights reserved.

Effortless Setup – Flexible Throughput – Precise Results

Quantify nucleic acids using outstanding data analysis algorithms, on a cycler with an easily interchangeable 96- or 384-well thermoblock.

- **Combine unsurpassed hardware with dedicated reagents and probes for pinpoint specificity and sensitivity.**
- **Discover and characterize mutations using powerful detection methods, including High Resolution Melting.**
- **Build custom panels or order complete single assays of your choice, using the free online RealTime ready Configurator.**

Experience the LightCycler® 480 Real-Time PCR System's ideal combination of accuracy, versatility, and speed.

Learn more at www.lightcycler480.com
and www.realtimeready.roche.com

Roche Diagnostics GmbH
Roche Applied Science
Werk Penzberg
82372 Penzberg, Germany




EXCEPTIONAL PERFORMANCE

Restriction Enzymes from New England Biolabs

With 35 years of experience in enzyme technology, New England Biolabs leads the industry in the study and improvement of restriction enzymes. NEB is pleased to introduce a line of High Fidelity (HF) enzymes that have been engineered for reduced star activity. HF enzymes bring a new level of flexibility when choosing reaction conditions, including volume, incubation time and buffer compatibility. Make NEB your first choice and experience the exceptional performance of this next generation in restriction enzyme technology.

Reduced star activity with HF Enzymes from NEB

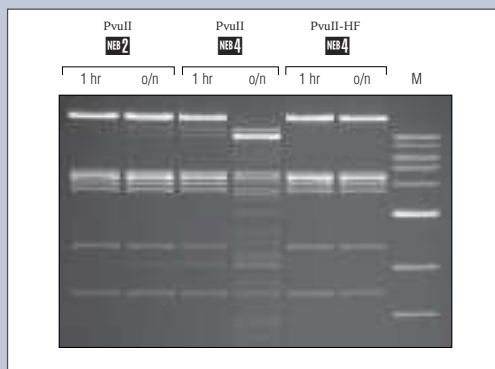

Recombinant PvuII, normally supplied in NEBuffer 2, exhibits star activity when used in extended incubations or in alternate buffers. Star activity is significantly reduced with PvuII-HF, even under extended incubation.

Advantages:

- **Selection** – More specificities than any other supplier
- **Convenience** – Optimal activity for over 160 enzymes in a single buffer **NEB 4**
- **Quality** – State-of-the-art production and stringent QCs
- **Innovation** – HF enzymes engineered for reduced star activity
- **Performance** – Guaranteed



CLONING & MAPPING

DNA AMPLIFICATION
& PCR

RNA ANALYSIS

PROTEIN EXPRESSION &
ANALYSIS

GENE EXPRESSION
& CELLULAR ANALYSIS

www.neb.com

Everyday needs. Everyday

answers.™



Easily find your favorite brands of reagents, kits, and benchtop devices at www.invitrogen.com/everyday.

Ambion®
RNA and RNA sample prep

Dynal®
Magnetic beads for separation

Gateway®
Cloning and protein expression

GIBCO®
Cell culture

Lipofectamine™
Transfection reagents

Molecular Probes®
Fluorescent dyes and probes

Novex®
Protein separation and blotting

SuperScript®
cDNA synthesis reagents

TaqMan®
Real-time PCR

TOP0®
PCR cloning

For research use only. Not intended for any animal or human therapeutic or diagnostic use, unless otherwise stated.

© 2009 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners. TaqMan is a registered trademark of Roche Molecular Systems, Inc. These products may be covered by one or more Limited Use Label Licenses (see the Invitrogen catalog or our website, www.invitrogen.com). By use of these products you accept the terms and conditions of all applicable Limited Use Label Licenses. A-085988 0609

Inspired Ägain

Who better to draw inspiration for the new ÄKTA™ avant system than from customers using the 30,000 ÄKTA systems already in use around the world? Well, you spoke and we listened. The new ÄKTA avant is faster — enabling quicker insights. It minimizes the chance of error, even while working at higher speeds. And it allows for direct, reliable scalability. At GE Healthcare, our focus is on helping scientists achieve even more, faster. It's a commitment we have in our genes. And all this is backed by the service, support, and investment in the future that being part of GE can bring.

Want to know more? Why not talk with us today. Visit www.gelifesciences.com/aktaavant

| ÄKTA™ | Amersham™ | Biacore™ | IN Cell Analyzer | Whatman™ | GE Service |

The New ÄKTA avant



imagination at work



It becomes you.

Introducing the 3500 Series Genetic Analyzer.

Get ready to make an amazing discovery: the new 8-capillary and 24-capillary 3500 Series Genetic Analyzers take DNA analysis to an entirely new level of performance. A level where your daily workflow seems like a natural extension of your own intuition. Where precise, quality-assured data inspires greater confidence. And where a new consumables design and intuitive software interface keep you current and in control.

Take a closer look, and you'll find the new 3500 and 3500xL Genetic Analyzers are like second nature. Which is our first priority when it's your data.

Discover the 3500 System at www.appliedbiosystems.com/3500Series



**Easy-to-Use
Consumables**



**Control at
Your Fingertips**



**Quality-Assured
Data**

AB applied
biosystems™

FOR RESEARCH USE ONLY. NOT FOR USE IN DIAGNOSTIC PROCEDURES.

© 2009 Life Technologies Corporation. All rights reserved. The trademarks mentioned herein are the property of Life Technologies Corporation or their respective owners.

For those who require IVD-marked devices, the 3500 Dx and the 3500xL Dx Genetic Analyzers and system accessories meet the requirements of the In Vitro Diagnostics Medical Devices Directive (98/79/EC). The 3500 Dx and 3500xL Dx systems are for distribution and use in specific European countries only. For more information about the 3500 Dx Series Systems, contact your Applied Biosystems representative.

WHERE THE BASIC SCIENCE

W
W



R
R
L
L
D
D

OF MEDICINE

Submit your work to *Science Translational Medicine* today!

This fall, AAAS and *Science* will launch *Science Translational Medicine*, a new journal focused on applications of basic research knowledge that will improve human health.

The goal of *Science Translational Medicine* is simple: to help the scientific community harness decades of progress in research at the basic level and translate these biological discoveries into medical advances. Take this opportunity to have your work recognized in this groundbreaking new journal.

Papers in the following areas will be reviewed and considered for publication:

- Animal & Human Studies
- Applied Physical Sciences
- Behavior
- Bioengineering
- Biomarkers
- Cancer
- Cardiovascular Disease
- Cell Culture
- Chemical Genomics/Drug Discovery
- Data Mining
- Drug Delivery
- Gene Therapy/Regenerative Medicine
- Imaging
- Immunology/Vaccines
- Infectious Diseases
- Medical Informatics
- Medical Nanotechnology
- Metabolism/Diabetes/Obesity
- Neuroscience/Neurology/Psychiatry
- Pharmacogenetics
- Policy
- Toxicology & Pharmacokinetics
- And other interdisciplinary approaches to medicine


To submit your work for consideration, please visit: www.submit2scitranslmed.org



www.ScienceTranslationalMedicine.org

INTEGRATING MEDICINE AND SCIENCE

Now you have a choice: 96 post-docs with Clark Electrodes. Or one Seahorse.

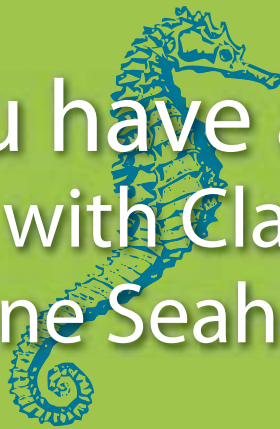

*Cellular
Bioenergetics
THE
Label-Free
Cell Based Assay*

Accelerate your cellular bioenergetics research without expanding your lab staff. All you need is one XF Extracellular Flux Analyzer from Seahorse Bioscience.

XF Extracellular Flux Analyzers from Seahorse Bioscience offer a new, better and faster way to profile cellular bioenergetics. Our breakthrough XF assays measure both mitochondrial respiration and glycolysis simultaneously, in real time--providing immediate insight into the effect of drug treatments, gene transfections and toxins on cell metabolism.

Our unique Transient Microchamber technology enables you to perform precise, label free, non-destructive measurements in 24 or 96 well microplates simultaneously. And with integrated drug delivery you can generate more physiologically relevant data in less time.

ATTEND OUR FREE WEBINAR SERIES ON CELLULAR BIOENERGETICS

Learn how leaders in metabolic and neurodegenerative disease, cancer, and in vitro toxicology are using our XF Extracellular Flux Analyzers to accelerate their research.



The new XF96 Extracellular Flux Analyzer

Register today at www.seahorsebio.com

 **Seahorse Bioscience™**

WHO NEEDS LABELS?

MACROMOLECULAR INTERACTION SANS LABELS

When it comes to monitoring biological interactions, whether of nucleic acids, proteins, small molecules, or cells, labels are essential for detection, right? Researchers certainly have embraced detection modalities based on tagging with radioisotopes, fluorescent dyes, and conjugated enzymes. Yet in many cases such molecular baggage is unnecessary, and possibly detrimental. But there is another way. Whether based on optical, mechanical, or electrical detection, label-free methods can probe molecular pas de deux in their native, unadulterated states. As a result, they are faster, simpler, and more physiological than their labeled counterparts. **By Jeffrey M. Perkel**

The research community is taking note. Says Guenter Gauglitz, whose lab does extensive work in this area, when asked why researchers should eschew tried-and-tested label-based methods for label-free: “I would pose the question the opposite way: Why bother with labeled methods?”

Tagging methods, says Gauglitz, of the **University of Tübingen**, Germany, have one very obvious disadvantage: “You have to label the systems.” That, in turn, can impact the “bioactivity” of the labeled molecule, for instance through misfolding, reduced mobility, or steric hindrance. Labeling reagents are also pricey, may tag different molecules with varying efficiency, and in some cases are incompatible with live cells.

Label-free approaches, in contrast, reduce biomolecular interactions to their most essential components: molecule A, molecule B, and a detection scheme to watch them in action, often in real time.

Researchers have, to some extent, caught on, but not in droves. “To date it’s been a niche technology,” says Matthew Cooper of the **University of Queensland**, Australia, editor of the recently published book *Label-Free Biosensors: Techniques and Applications* (Cambridge University Press, 2009). “It’s been a classic early-adopter play.”

Pharmaceutical companies have been more enthusiastic, pressing label-free systems into service for fragment-based screening, lead optimization, and affinity ranking, among other applications.

At biopharmaceutical firm **Amgen**, as head of the company’s protein-protein interaction group, principal scientist Ching Chen uses three different label-free platforms. Those systems, from **GE Healthcare**, **ForteBio**, and **Sapidyne Instruments**, are used, she says, “in a complementary way,” for such applications as primary screening, epitope mapping, and affinity determination.

Label-free technologies provide “a much more straightforward answer [than labeled approaches],” she says. “They provide much better confidence when interpreting the data and give us confidence and speed in the selection of molecules.”

Three Essentials

All label-free (and for that matter, labeled) technologies, regardless of their operating principle and configuration, have at their heart three fundamental components, says Cooper: a receptor, a transducer, and an analyzer.

The receptor, of course, is what makes a biosensor, a biosensor—it specifically captures the molecule, cell, or pathogen of interest from a biological fluid or sample. The transducer, says Cooper, is the “detection modality, what converts that biological event into a photometric, electrochemical, electrical, or acoustical change.” That change, which is proportional to the extent of binding, is then quantified by the analyzer.

Observe how these pieces work together in one of the most popular label-free approaches, surface plasmon resonance (SPR). [continued >](#)



“Label-free technologies provide “a much more straightforward answer” [than labeled approaches].”


Look for these Upcoming Articles

Epigenetics — October 9

RNAi — October 16

Inclusion of companies in this article does not indicate endorsement by either AAAS or Science, nor is it meant to imply that their products or services are superior to those of other companies.

Proteomics

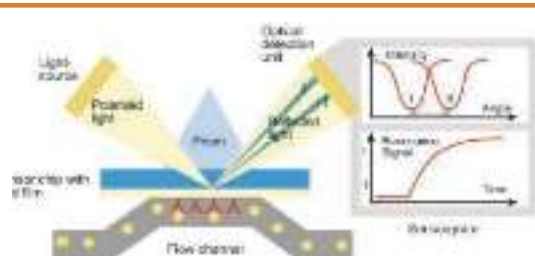

As implemented by GE Healthcare, an SPR biosensor comprises a thin gold film (about 50 nm thick) mounted on the flat face of a prism. On the side of the gold film facing away from the prism is a surface on which capture molecules (such as antibodies or DNA) are immobilized. This surface (the transducer) is exposed to a fluidics interface, such that reagents—whether capture molecules, analytes, or buffer—can flow over and interact with the biosensor surface (see figure above).

Within the instrument itself, the biosensor, light source, and detector form the vertices of a triangle, an arrangement called the Kretschmann configuration. As polarized light passes through the prism at an angle, it strikes the gold and reflects back through the prism at the “specular” angle—that is, as if the gold were a mirror. For the most part, the gold really is a mirror, and over a wide range of wavelengths (or angles), the light hitting the detector (a photodiode that measures light intensity) is 100 percent of that which left the source. But at one particular combination of wavelength and angle, reflection drops nearly to zero—a phenomenon called the “SPR dip” or “SPR resonance.”

That’s because under those precise conditions, “the light resonantly couples with the conduction electrons near the surface of the gold film to create clouds of oscillating electrons called a surface plasmon,” says Charlie Campbell, the West Professor of Chemistry at the **University of Washington**. “When in resonance, the momentum of this plasmon parallel to the surface matches the parallel momentum of the incoming light; since the latter is a function of the refractive index of the material under analysis in contact with the gold surface, the angle [at fixed wavelength] or wavelength [at fixed angle] of the dip can be used to monitor the refractive index of that material.”

In practice, this index of refraction varies ever so slightly as molecules bind to the sensor surface. So as intermolecular binding events unfold, the wavelength at which the SPR “dip” occurs will change, and the magnitude of that change is a function of the extent of binding.

The “refractive index change is directly proportional to the mass change,” says Stefan Lofas, principal scientist for the **Biacore** line at GE Healthcare Life Sciences. “So actually, you can say that we quantitatively measure the change in mass when something binds to the surface.”

And, because SPR devices are often coupled to flow cells that can deliver analytes or buffer, they can measure those events in real-time (in contrast to systems such as ELISA, which collect end-point measurements). That makes SPR systems particularly useful for collecting kinetics data—and explains their popularity in the pharmaceutical industry as well.

Biopharmaceutical firm **Xencor** runs its Biacore 3000 “24/7” in its antibody development pipeline, says John Desjarlais, vice president of research; **MedImmune’s** SPR instruments (from Biacore and **Bio-Rad**) have been used for affinity kinetic characterization, epitope competition studies, and characterizing concurrent binding of

multispecific antibodies, says Herren Wu, vice president of research and development.

“In the past 10 years, almost all compounds have gone through the Biacore measurement,” he says.

By all accounts, Biacore has been a dominant force in the label-free market in general and SPR in particular for nearly two decades. “Last fall, we celebrated the 10,000th publication with Biacore-related data,” Lofas says. Yet competition is growing. **Bio-Rad’s** ProteOn system can measure 36 interactions simultaneously via a criss-crossing array of 6x6 fluidic channels, while **Fujifilm Life Science’s** AP-3000 can probe six targets at once. An in-development system from **Plexera**, on the other hand, applies an SPR variant called imaging SPR (or SPR microscopy) to monitor interactions on arrays of up to several thousand spots by imaging light intensity using CCD cameras.

Yet all these techniques suffer from the same basic shortcomings, namely, that they are flow-based and mounted on a prism. As a result, SPR is generally incompatible with standard microtiter-plate formats, including much of the high throughput drug discovery workflow.

Optical Alternatives

Several companies have overcome this problem by devising SPR alternatives in standard microplate labware. Though flow-based systems have the advantage of providing kinetics data, microplate-based systems can accommodate a wider variety of sample types, from cruder samples that would gum up microfluidic channels, to live cells that are currently too large for them.

The majority of these options—like SPR—rely on evanescent waves. As described in a review by Cooper in *Nature Reviews Drug Discovery* (1:515, 2002), the “evanescent wave phenomenon” occurs when “total internal reflection of light at a surface–solution interface produces an electromagnetic field, or evanescent wave, that extends a short distance (~100–200 nm) into the solution.” That wave effectively probes the interface—and is in turn “tuned” by it in a dose-dependent fashion—allowing detection.

Just how that detection occurs depends on what aspect of the outbound light the detector looks at. **Corning’s** Epic system probes the index of refraction on its biosensor surface using a “resonant waveguide grating structure,” says Mark Krol, technology manager at Corning.

“From a phenomenological standpoint, [Epic] is very similar to SPR, but the optical physics are very different,” Krol explains.

As with SPR, molecular interaction events on the sensor surface change its index of refraction. But rather than monitoring the SPR dip, Epic scans for a concomitant change in reflected wavelength, whose magnitude is a function of the degree of binding.

SRU Biosystems’ BIND performs a similar measurement using an optical device called a photonic crystal.

“A photonic crystal is just a fancy way of saying we have a material with two different refractive indexes in it, a high refractive index and a low refractive index, that alternate from high to low to high to low in a periodic fashion,” explains company co-founder Brian Cunningham of the **University of Illinois, Urbana-Champaign**.

In this case, the crystal is illuminated with near-infrared light, almost all of which passes through the sensor. But one wavelength resonates with the crystal and is reflected back. “At that one wavelength, the resonant wavelength, the photonic crystal structure acts as a perfect mirror,” Cunningham explains. “So it reflects very efficiently at that one wavelength, and all other wavelengths pass through.” When a biomolecular interaction occurs on the sensor surface, that wavelength shifts, producing a measurable readout. **continued >**



Lead the next generation of pharmaceutical science.

Discover the Answers that Matter.

Eli Lilly and Company is a leading, innovation-driven pharmaceutical corporation with a firm commitment to help people live longer, healthier and more active lives by making breakthroughs in medicines and treatments. To find the next generation of drugs, we know we need to use next generation technologies and employ the best of the biopharmaceutical industry.

Established in 2002, Lilly Singapore Centre for Drug Discovery (LSCDD) is now expanding its capabilities to discover and develop new medicines more productively, in the areas of cancer and metabolic disorders. We form a network of drug development partners, and through innovative data integration approaches, discover and apply biomarker and patient-tailoring solutions.

Located in the exciting Singapore Biopolis, LSCDD's multidisciplinary and multicultural team is working to redefine the leading edge. We are looking for outstanding individuals with demonstrated industry experience to fill the following positions:

- **Group Leader (Cancer Biology)**
- **Group Leader (Diabetes Biology)**
- **Bioinformatics Group Leader (Systems Biology)**
- **Epigenetics Platform Biology Leader**
- **Bioinformatics Manager (Integrative Computational Sciences)**
- **Principal Statistician (Integrative Computational Sciences)**
- **Senior Statistical Geneticist (Integrative Computational Sciences)**
- **Senior Medicinal Chemist**
- **Postdoctoral Fellow, Bioinformatics (Integrative Computational Sciences)**
- **Postdoctoral Fellow, Kinase (Oncology)**
- **Postdoctoral Fellow (Systems Biology)**

Please log on to www.lscdd.lilly.com.sg to find out more about these positions.

www.lscdd.lilly.com.sg

Lilly

Answers That Matter.



Life Even Better™

Life Technologies is providing tools for quick, accurate diagnosis during public health emergencies. In the 2009 H1N1 outbreak, we hastened detection of influenza A by providing instrument systems and components to public health labs across the U.S., enabling urgent CDC-consistent virus testing. This is just one example of how we enable people to accelerate scientific exploration, driving to discoveries and developments that make life even better.

We believe in the power of science to transform life and the power of people to make it happen. Join us in discovering the infinite surprises and solutions that science—and you—can reveal.

www.lifetechnologies.com/careers
EOE

life
technologies™

AB applied biosystems™ | **i** invitrogen™

Cunningham compares the effect to a guitar string. “If you strum it, you get a certain wavelength out in terms of the sound it makes,” he says. “If you paint the string and strum it again, you will get a different resonance, a different sound. In the same way, when mass is added to this sensor, it reflects the light at a different wavelength.”

Maven Biotechnologies’ in-development Polaron instrument measures changes in light polarization, says chief scientific officer Shane Dultz. Based on a technology called LFIRE (label-free internal reflection ellipsometry), Polaron’s is an imaging approach that measures changes in the elliptical polarization of light as it interacts with the sensor surface to chronicle biomolecular interactions in miniaturized protein microarrays.

“It turns out that the phase shift, or the degree to which the light’s elliptical polarization is changing, is proportional to the number of molecules per square micron on that surface,” Dultz says.

Though still in development, Polaron will image arrays in the bottom of 96-well glass bottom plates, Dultz says, with up to 700 protein spots per well. Each spot, he adds, can be analyzed either continuously or on an endpoint basis.

Silicon Kinetics (SKi Pro system), **BiOptix**, and **ForteBio** (Octet) employ an alternative strategy called interferometry. Used, among other things, to calculate the distance to astronomic objects, interferometry may also be used to measure the thickness of the molecular layer on the biosensor surface as biochemical events occur.

In ForteBio’s case, the biosensor is a narrow tube (or fiber tip) with two parallel surfaces perpendicular to its long axis: an invariant reference surface and a sensor surface, to which capture molecules (such as antibodies) are bound. When light is directed down the barrel of the tube, it bounces off both surfaces on its way back to

the detector, producing a composite signal that is a reflection of both waves.

“The two reflections interfere constructively and destructively and there is a net signal seen at the detector,” explains Sriram Kumaraswamy, product manager of ForteBio.

As molecules bind to the biosensor surface, that net signal pattern changes, producing a difference in the resulting spectrum whose magnitude depends on the size and number of molecules bound.

Electrical and Mechanical Biosensing

Though most commercial systems employ optical detection, there are other ways to detect interactions sans labels. As Cunningham puts it, “All optical biosensors take advantage of the fact that biomolecules all have a dielectric coefficient or refractive index that is higher than that of water. But what other properties does a molecule have?”

One property, of course, is mass, and several mechanical approaches exist to probe that, says Rashid Bashir, the Bliss Professor of Electrical and Computer Engineering and Bioengineering at the University of Illinois at Urbana-Champaign. One such approach is “resonance mass sensing.”

Many structures—quartz crystals, for instance—have a characteristic vibrational frequency, Bashir explains. Because that frequency depends on mass, if the device is small enough, “as I add a mass to it, I can detect that added mass by the change in the resonant frequency of the resonator.”

TTP LabTech’s RAPid 4 system is a quartz crystal microbalance device based on this principle.

An alternative strategy, stress sensing, uses mechanical stress to follow biochemical events on the surface of, for instance, a cantilever (a kind of nanoscale diving board). “If you perform a biochemical reaction on only one surface of these cantilevers, then that change in the local surface energy will result in a change in surface stress, and the cantilever will actually bend,” Bashir explains. As with other label-free methods, the degree of that flexion depends on the extent of macromolecular interaction, and can be quantified by bouncing a laser off the cantilever’s tip.

Others have devised label-free sensing methods based on electronic fluctuations resulting from macromolecular interactions. According to Bashir, systems have been devised based on impedance or electrochemical sensing using electrodes in fluid or potentiometric sensing using field-effect transistors or silicon nanowires. **Acea Biosciences** and **MDS Analytical Technologies** have commercialized electrical impedance-based systems, which measure changes in current flow between electrodes when cells deposited between the electrodes change shape, adherence, or mass distribution as a result of, for instance, G protein-coupled receptor stimulation.

Ultimately each strategy, whether optical, mechanical, or electrical, has its pros and cons. SPR, for instance, has decades of user know-how behind it. Electrical approaches may be most amenable to miniaturization, says Bashir, as they require neither optics nor vibrational detectors, whereas electrical and mechanical strategies have the advantage of being immune to optical distortions at the biosensor surface.

Whichever molecular property you focus on, once you get over the learning curve, you might find that once you go label-free, there’s no going back.

Jeffrey M. Perkel is a freelance science writer based in Pocatello, Idaho.

DOI: 10.1126/science.opms.p0900037

Featured Participants

Acea Biosciences
www.aceabio.com

Amgen
www.amgen.com

Biacore
www.biacore.com

BiOptix
www.biopix.com

Bio-Rad
www.bio-rad.com

Corning
www.corning.com/lifesciences

ForteBio
www.fortebio.com

GE Healthcare
www.gehealthcare.com

Maven Biotechnologies
www.mavenbiotech.com

MDS Analytical Technologies
www.moleculardevices.com

MedImmune
www.medimmune.com

Sapidyne Instruments
www.sapidyne.com

Silicon Kinetics
www.siliconkinetics.com

SRU Biosystems
www.srubiosystems.com

TTP Labtech
www.ttplabtech.com

University of Illinois, Urbana-Champaign
illinois.edu

University of Queensland
www.uq.edu.au

University of Tübingen
http://www.uni-tuebingen.de/en/landingpage.html

University of Washington
www.washington.edu

Xencor
www.xencor.com

New Products

**Scanning Microplate Reader with Spectrophotometer**

The compact Zenyth 200 is a high-performance scanning microplate reader with the added benefit of a built-in ultraviolet/visible light spectrophotometer. It is designed to meet the needs of life science laboratories that need a fast, accurate plate reader with the flexibility to handle multiple plate formats from cell culture plates (6-well to 48-well) to 96-well and 384-well plates. The Zenyth 200's advanced optics are capable of endpoint, kinetic, multichromatic, spectral, and position scanning, covering the range from 190 nm to 1,000 nm. The temperature can be controlled for kinetic studies. It is available in two versions: a personal-computer-controlled version and a stand-alone version with built-in routines for nucleic acid quantitation and characterization and integrated Windows CE for easy data transfer.

Biochrom

For information +44-(0)-1223-423723

www.biochrom.co.uk

Workstation for Transporter Research

The Surfe2r Workstation 5 is a semiautomated, benchtop system designed for basic research and assay development that features proprietary Surface Electrogenic Event Reader technology. The cell-free, biosensor-based system enables direct, precise measurement of transporter activity free of molecular labels, providing in-depth analysis of transporter function. The Workstation 5 provides one sensor position for a single sensor and up to 84 sample positions for buffers. The Surface Electrogenic Event Reader technology does not make use of entire cells but instead measures the charge displacement from membrane preparations that are isolated from cell culture and other sources. Membranes are coupled, via an alkaline thiol and lipid bilayer, to gold-coated electrodes that detect changes in the electrical signal, without the use of reporter molecules.

IonGate Biosciences

For information +49-69-305-85400

www.iongate.net

Label-Free Interaction Analysis

The RAPid 4 is a flow-based system that performs label-free interaction analysis through the use of novel Resonant Acoustic Profiling technology. The detection of molecular interactions through acoustic technology offers advantages over optical biosensors. It reduces the need to purify samples and generates accurate kinetic, affinity, and concentration measurement from complex biological mixtures, such as blood, serum, cell culture supernatants, and periplasmic extracts. The RAPid 4 is suitable for use in life science and pharmaceutical laboratories for biotherapeutics development and protein expression analysis. It exploits the piezoelectric effects of quartz crystal resonators for label-free detection of molecular interactions in real time.

TTP LabTech

For information 617-494-9794

www.ttplabtech.com

Label-Free Instruments

Two next generation, high throughput, label-free instruments have been launched. The Octet RED384 is for protein, peptide, small molecule, and fragment screening. The Octet QK384 is for protein and antibody assays. Both instruments allow 384-well detection, 16-channel simultaneous readout, biosensor regeneration and reracking, and automation capabilities for biotherapeutic and pharmaceutical drug discovery assays. The new instruments increase the cost efficiencies of running ForteBio's la-

bel-free Dip and Read biosensor assays in high throughput mode, while speeding the time to results for kinetic characterization and quantitation measurements. The systems support two 384-well or 96-well plate positions. A full complement of samples can be run in one plate, while the second plate holds buffers, diluents, and controls.

ForteBio

For information 650-289-6800

www.fortebio.com

Cell Migration Assay

The Oris Cell Migration Assay allows for real-time imaging and can be used on high throughput and high content imaging instruments. Unlike traditional cell migration assays using a Boyden chamber device, there is no transmembrane or cell culture insert obstructing live images of the cells as they migrate. The assay features a specially designed, cell-seeding stopper inside each microplate well that restricts cells from seeding into a central detection zone as cells are being dispensed. After cell attachment has occurred, the stoppers are removed, allowing cells to migrate into the detection zone for analysis using a fluorescence plate reader, microscope, or high content imaging system. The assays are available with tissue culture treated plates, collagen I coated plates, fibronectin coated plates, or a tricoated version with all three surfaces.

AMS Biotechnology

For information +44-1235-828200

www.amsbio.com

Protein and Fragment Assays

Two new label-free assay modules—for high-sensitivity protein quantitation and fragment screening—are available for use on the Octet instrument platform. The new modules are available via software upgrades. The high-sensitivity protein assay module can detect low nanogram per milliliter levels of analytes in quantitation mode. It can replace fluorescence and luminescence enzyme-linked immunosorbent assays to provide results within an hour, enabling follow-up research to be conducted the same day. The fragment screening assay extends the small molecule analysis capabilities available on the Octet RED system by enabling kinetic characterization and screening of fragment libraries that include compounds in the 150 to 300 dalton range.

ForteBio

For information 650-289-6800

www.fortebio.com

Electronically submit your new product description or product literature information! Go to www.sciencemag.org/products/newproducts.dtl for more information.

Newly offered instrumentation, apparatus, and laboratory materials of interest to researchers in all disciplines in academic, industrial, and governmental organizations are featured in this space. Emphasis is given to purpose, chief characteristics, and availability of products and materials. Endorsement by *Science* or AAAS of any products or materials mentioned is not implied. Additional information may be obtained from the manufacturer or supplier.

Do it Right the *First Time* PCR that Works!

Takara offers a wide selection of PCR enzymes for most popular applications.

• High Fidelity PCR

PrimeSTAR[®] HS DNA Polymerase for 10X higher fidelity than *Taq* as well as extended lengths up to 20 kb.

• Real Time PCR

SYBR[®] Premix Ex Taq[™] (Perfect Real Time) for sensitive qPCR detection using SYBR[®] Green I.

SYBR[®] Premix Ex Taq[™] II (Perfect Real Time) for improved reaction specificity and performance.

Premix Ex Taq[™] (Perfect Real Time) for specific qPCR using probe detection.

All three premixes are compatible with a variety of real-time PCR instruments.

• Multiplex PCR

Taq Hot Start DNA Polymerase for excellent sensitivity and reproducibility in multiplex PCR.

LA Taq[™] HS DNA Polymerase for challenging multiplex PCR.

• High-Speed PCR

SpeedSTAR[™] Hot Start DNA Polymerase for high speed, high sensitivity PCR.

• Hot-Start PCR


Ex Taq[™] HS DNA Polymerase for increased sensitivity and high yield.

Taq HS DNA Polymerase for increased specificity and reduced background.

LA Taq[™] HS DNA Polymerase for long PCR.

• Long PCR

LA Taq[™] DNA Polymerase for superior amplification of long and GC-rich templates.

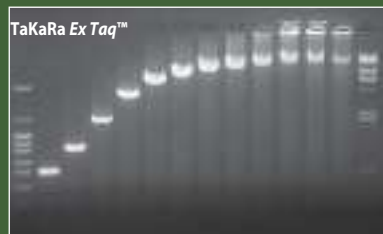

Takara's robust, sensitive enzymes yield superior results on templates of varying sizes. *Ex Taq[™]* produces excellent amplification of fragments from 0.5 kb to 20 kb.

PrimeSTAR[™], SpeedSTAR[™], TaKaRa Ex Taq[™] and TaKaRa LA Taq[™] are trademarks of Takara Bio Inc. LA PCR technology is covered by U.S. Patent No. 5,436,149 issued to Takara Bio Inc. Purchase of this product includes an immunity from suit under patents specified in the product insert to use only the amount purchased for the purchaser's own internal research. No other patent rights (such as 5' Nuclease Process patent rights) are conveyed expressly, by implication, or by estoppel. Further information on purchasing licenses may be obtained by contacting the Director of Licensing, Applied Biosystems, 850 Lincoln Centre Drive, Foster City, California 94404, USA. Takara Bio's Hot-Start PCR related products are licensed under U.S. patent 5,338,671 and 5,587,287 and corresponding patents in other countries. SYBR[®] is a registered trademark Molecular Probes.

TaKaRa

For more information
www.takara-bio.com

Japan:
Takara Bio Inc.
+81 77 543 7247
www.takara-bio.com

USA:
Clontech Laboratories
A Takara Bio Company
888-251-6618
www.takara-bio.us

Europe:
Takara Bio Europe S.A.S.
+33 1 3904 6880
www.takara-bio.eu

China:
Takara Biotechnology
(Dalian) Co., Ltd.
+86 411 8764 1681
www.takara.com.cn

Korea:
Takara Korea
Biomedical Inc.
+82 2 2081 2525
www.takara.co.kr

amplifamily.

A family of solutions for better PCR.

You talked PCR. We listened. Our 20-year collaboration with you has led to the development of a family of amplification products that delivers the best combination of performance and value.

Thanks to you, our comprehensive portfolio of thermal cyclers, qPCR systems, software, reagents, plates, and tubes builds on our strong lineage and continues to define what you expect in PCR innovations.

Reap the benefits of the solutions you've helped us create. See what rethought PCR can do for you.

For more information on our complete suite of products and to request a sample of our newest qPCR reagent, please visit us at www.bio-rad.com/ad/amplifamily/.

Research. Together.

